

Two-color fluorescent in situ hybridization using chromogenic substrates in zebrafish

Jennifer A. Schumacher, Emma J. Zhao, Matthew J. Kofron, and Saulius Sumanas

Division of Developmental Biology, Cincinnati Children's Hospital Medical Center, Cincinnati, OH

BioTechniques 57:254-256 (November 2014) doi 10.2144/000114229

Keywords: fluorescent in situ hybridization; FISH; two color; NBT; BCIP; Vector Red; zebrafish; confocal

Two-color fluorescent in situ hybridization (FISH) is a widely used technique for comparing relative gene expression patterns. Current two-color FISH protocols are not ideal for detecting weakly expressed transcripts or monitoring signal strength and background levels during the course of the reaction. Here we describe an improved FISH protocol using the conventional highly sensitive chromogenic substrates nitro blue tetrazolium (NBT)/5-bromo-4-chloro-3-indolyl phosphate (BCIP) and Vector Red in zebrafish embryos. This protocol substantially improves on existing FISH techniques by combining the advantages of long reactivity of alkaline phosphatase, chromogenic monitoring of both developing reactions, and the ability to perform subsequent high-resolution fluorescent imaging. Although tested in zebrafish, a similar approach is expected to be applicable to ISH in any model organism.

In situ hybridization (ISH) is the method of choice for visualizing the distribution of transcripts in developing embryos. Simultaneous detection of two transcripts is commonly used to determine overlap in expression domains. Two-color ISH protocols have been reported in various organisms using simultaneous hybridization with digoxigenin (DIG)- and fluorescein (FL)-labeled probes followed by sequential visualization with alkaline phosphatase (AP) chromogenic substrates, including NBT/BCIP and Fast Red (1,2). The major advantage of chromogenic reactions is the ability to monitor the AP reaction to control signal strength and background levels. Furthermore, the AP reaction has a long productivity time to

help detect weakly expressed transcripts. NBT/BCIP, which produces a blue-purple precipitate, is generally the substrate of choice for chromogenic ISH because of the strong signal and low background levels the reaction generates. However, with two-color ISH, the darker NBT/BCIP substrate often masks the lighter Fast Red substrate, making it difficult to determine if transcripts are co-expressed in the same cell.

To address this problem, multiplex FISH techniques that rely on horseradish peroxidase (POD) detection followed by sequential tyramide signal amplification (TSA) visualization have been developed (3–5). Although these protocols allow for high-resolution imaging to examine

overlapping expression in single cells, they do not allow for monitoring of the POD-TSA developing steps. The tyramide substrate is fluorescent on its own, and extensive washing is required to visualize signal. Additionally, because the enzymatic reaction lasts for only a short time, detecting weakly expressed transcripts can be challenging. A two-color FISH protocol for combined AP-Fast Blue and POD-TSA was recently described; however, this protocol still requires one probe to be visualized without monitoring the developing step (6).

We report here a two-color FISH protocol in zebrafish embryos that takes advantage of the fluorescent properties of the NBT/BCIP and Vector Red substrates. These substrates were chosen based on non-overlapping emission wavelengths; NBT/BCIP fluoresces in the near-infrared range (7), and Vector Red is visible with Texas Red or rhodamine filter sets, similar to Fast Red (6). Importantly, this protocol combines the advantages of long AP reactivity, chromogenic monitoring of both developing reactions, and the ability to perform subsequent high-resolution fluorescent imaging.

Twenty somite stage zebrafish embryos were fixed in 4% paraformaldehyde (PFA) following standard procedures (8). They were dehydrated through an ethanol (EtOH) series and stored at -20°C . On day 1, embryos were rehydrated, washed 3 times in PBT [1 \times PBS (13.7mM NaCl, 0.27 mM KCl, 0.43 mM Na_2HPO_4 , 0.14 mM KH_2PO_4) with 0.2% bovine serum albumin, 0.2% Tween 20], and then incubated in prehybridization buffer [50% formamide, 5 \times SSC (for 20 \times SSC stock: 3M NaCl, 0.3 M sodium citrate, pH to 7.0), 50 $\mu\text{g}/\text{mL}$ heparin, 5 mM EDTA, 0.5 mg/mL torula yeast RNA, 0.1% Tween 20, pH 6.0 with citric acid] for 2 h at 65°C . Embryos were incubated in probes diluted in prehybridization buffer at 65°C overnight. Two probes were applied simultaneously; ideally the stronger probe was labeled with FL and the weaker probe with DIG. We tested the combination of *myeloperoxidase* (*mpx*)-DIG, a relatively weak probe that labels neutrophilic granulocytes

METHOD SUMMARY

We describe an improved fluorescent in situ hybridization (FISH) protocol using the conventional highly sensitive chromogenic substrates nitro blue tetrazolium (NBT)/5-bromo-4-chloro-3-indolyl phosphate (BCIP) and Vector Red in zebrafish embryos. This protocol substantially improves on existing FISH techniques by combining the advantages of long reactivity of alkaline phosphatase, chromogenic monitoring of both developing reactions, and the ability to perform subsequent high-resolution fluorescent imaging.

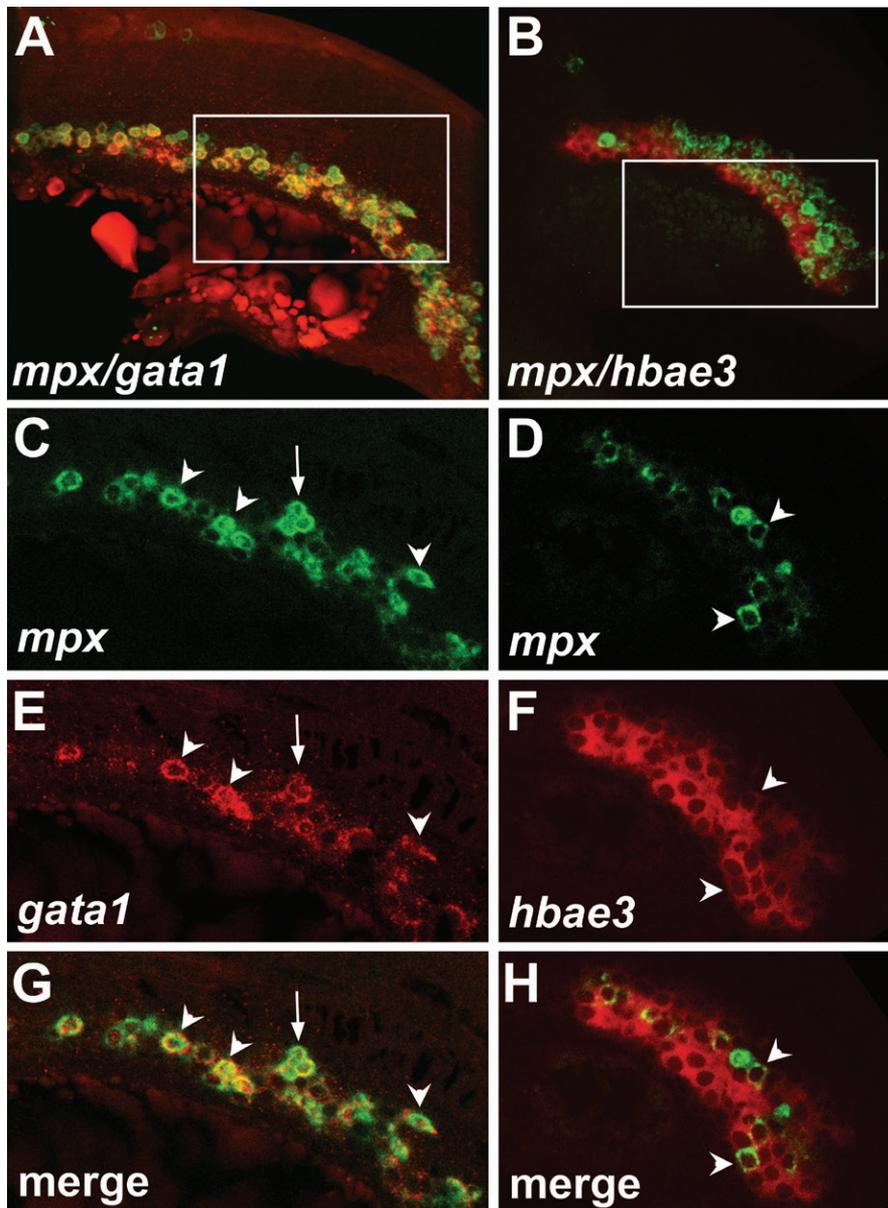


Figure 1. Fluorescent imaging of NBT/BCIP and Vector Red in situ hybridization (ISH). (A) Volume reconstruction of *mpx* (green) and *gata1* (red) expression. (B) Volume reconstruction of *mpx* (green) and *hbae3* (red) expression. (C,D) *mpx* transcript localization in the inner cell mass (ICM) visualized by nitro blue tetrazolium (NBT)/5-bromo-4-chloro-3-indolyl phosphate (BCIP). (E) *gata1* transcript localization in erythrocyte precursors visualized by Vector Red. (F) *hbae3* transcript localization in erythrocytes visualized by Vector Red. (G) Merged channels show overlapping expression of *mpx* and *gata1* in the majority of cells. (H) Merged channels show overlapping expression of *mpx* and *hbae3* in a subset of *hbae3*-positive cells. Volume reconstructions were performed using Imaris software (Bitplane, Zurich, Switzerland). Images in (C-H) are single confocal slices. Arrowheads indicate cells with expression of both transcripts. Arrows indicate cells with lower expression of *gata1*. All embryos are at the 20 somite stage, anterior to the left.

and exhibits variable expression levels in the inner cell mass (ICM) of zebrafish embryos (9) and either *gata1*-FL or α embryonic hemoglobin (*hbae3*)-FL, both of which are expressed in erythrocytes in the ICM (10,11). Both anti-DIG-AP and anti-FL-AP (Roche Diagnostics, Indianapolis, IN) were preabsorbed with acetone powder as previously described (8).

On day 2, embryos were washed with the following series of buffer washes at 65°C: 75%, 50%, and 25% prehybridization buffer diluted in 2 \times SSC for 15 min each, 2 \times SSC for 15 min, and then twice in 0.2 \times SSC for 30 min each. Embryos were washed in a dilution series of 0.2 \times SSC:PBT (3:1, 1:1, 1:3, PBT) for 5 min each at room temperature. During the washes, the 10 \times anti-DIG-

BLOW UP THE BARRIERS TO YOUR NEXT-GEN SEQUENCING



Next-gen sample QC is now hassle free.

FULLY AUTOMATED FRAGMENT ANALYZER™ DOES IT ALL.

- Assesses quality and quantity (size and concentration)
- Resolves fragments from 25 bp to 5,000 bp
- Sizes fragments up to 20,000 bp for PacBio sequencers
- Also analyzes gDNA and RNA

No chips. No tapes. No compromises.

More at AATI-US.COM



AP was spun down and the supernatant was diluted to 1× in 2% lamb serum in PBT. Embryos were incubated overnight at 4°C in anti-DIG-AP. We performed the DIG-AP reaction first because the first AP reaction has higher sensitivity and provides more sensitive detection of the weaker probe.

On day 3, embryos were washed 6 times in PBT, once in AP buffer (100 mM Tris pH 9.5, 100 mM NaCl, 50 mM MgCl₂, 0.1% Tween 20), and developed in 400 μL of developing buffer (4.5 μL of 50 mg/mL NBT and 3.5 μL of 50 mg/mL BCIP per 1 mL of AP buffer). When the desired signal intensity was reached, the reaction was stopped with several PBT washes. The anti-DIG-AP was inactivated by fixing in 4% PFA for 1 h at room temperature. Embryos were similarly processed for anti-FL-AP, except that AP buffer was replaced with 0.2 M of Tris pH 8.5 with 0.1% Tween 20. The FL probe was developed using Vector Red substrate (Vector Labs, Burlingame, CA) according to the manufacturer's instructions.

Dehydrating embryos in EtOH overnight after developing markedly reduced background during fluorescent imaging, particularly for the NBT/BCIP substrate. Imaging is ideally performed within a few days after processing for the highest signal-to-noise ratio. Embryos were mounted in 0.6% low melting point agarose overlaid with embryo media. Confocal images were acquired using a Nikon A1R si STORM inverted microscope equipped with an Apo LWD 40× /1.15 NA microscope objective (Nikon Instruments Inc., Melville, NY). NBT/BCIP fluorescence was excited using a 647-nm diode laser set to 100 mW intensity, detected with a 740-nm long pass filter, and then imaged with a high-quantum efficiency (QE) bi-alkali photomultiplier tube with 10% QE at 800 nm, similar to a previously reported protocol (7). Vector Red fluorescence was imaged with an excitation

wavelength of 561 nm and detected with a 595/50 nm emission filter.

Volume reconstructions of confocal Z-stacks show that the expression domain of *mpx* clearly overlaps with both *gata1* and *hbae3* (Figure 1, A and B). In confocal sections, individual cells expressing *mpx*, *gata1*, or *hbae3* could be clearly identified in the ICM region, with very little background fluorescence (Figure 1C–F). There was definitive overlap of *mpx* and *gata1* expression in the majority of *mpx*-expressing cells (Figure 1, C, E, and G) as previously reported (12). Some cells had strong expression of *mpx* and weaker expression of *gata1*, but with the ability to monitor both developing reactions, followed by fluorescent imaging, clear overlap can be seen. Only a subset of *hbae3*-positive cells also express *mpx* (Figure 1, D, F, and H), indicating Vector Red is not being detected by the NBT/BCIP filter.

In conclusion, this FISH protocol allows researchers to examine the distribution and overlap of two gene transcripts using the highly sensitive AP substrates NBT/BCIP and Vector Red. The ability to monitor both developing reactions and the higher sensitivities of the NBT/BCIP and Vector Red substrates will substantially improve the ability to visualize weak or variably expressed genes.

Author contributions

J.A.S. and E.J.Z. performed experiments; J.A.S., E.J.Z., and S.S. analyzed data; M.J.K. assisted with confocal imaging; and J.A.S. and S.S. prepared the manuscript.

Acknowledgments

Research was supported by Cincinnati Children's Research Foundation and NIH R01 HL107369 award to S.S. J.A. Schumacher was supported by NIH T32HL00738 fellowship, and E.J. Zhao was supported

by AHA-Great Rivers Affiliate Summer Undergraduate Research Fellowship 12UFEL9990000. This paper is subject to the NIH Public Access Policy.

Competing interests

The authors declare no competing interests.


References

1. Jowett, T. and L. Lettice. 1994. Whole-mount in situ hybridizations on zebrafish embryos using a mixture of digoxigenin- and fluorescein-labelled probes. *Trends Genet.* 10:73-74.
2. Hauptmann, G. and T. Gerster. 1994. Two-color whole-mount in situ hybridization to vertebrate and Drosophila embryos. *Trends Genet.* 10:266.
3. Brend, T. and S.A. Holley. 2009. Zebrafish whole mount high-resolution double fluorescent in situ hybridization. *J. Vis. Exp.* 25:e1229.
4. Clay, H. and L. Ramakrishnan. 2005. Multiplex fluorescent in situ hybridization in zebrafish embryos using tyramide signal amplification. *Zebrafish* 2:105-111.
5. Kosman, D., C.M. Mizutani, D. Lemons, W.G. Cox, W. McGinnis, and E. Bier. 2004. Multiplex detection of RNA expression in Drosophila embryos. *Science* 305:846.
6. Lauter, G., I. Soll, and G. Hauptmann. 2011. Two-color fluorescent in situ hybridization in the embryonic zebrafish brain using differential detection systems. *BMC Dev. Biol.* 11:43.
7. Trinh le, A., M.D. McCutchen, M. Bonner-Fraser, S.E. Fraser, L.A. Bumm, and D.W. McCauley. 2007. Fluorescent in situ hybridization employing the conventional NBT/BCIP chromogenic stain. *Biotechniques* 42:756-759.
8. Jowett, T. 2001. Double in situ hybridization techniques in zebrafish. *Methods* 23:345-358.
9. Bennett, C.M., J.P. Kanki, J. Rhodes, T.X. Liu, B.H. Paw, M.W. Kieran, D.M. Langenau, A. Delahaye-Brown, et al. 2001. Myelopoiesis in the zebrafish, *Danio rerio*. *Blood* 98:643-651.
10. Detrich, H.W., 3rd, M.W. Kieran, F.Y. Chan, L.M. Barone, K. Yee, J.A. Rundstadler, S. Pratt, D. Ransom, and L.I. Zon. 1995. Intraembryonic hematopoietic cell migration during vertebrate development. *Proc. Natl. Acad. Sci. USA* 92:10713-10717.
11. Brownlie, A., C. Hersey, A.C. Oates, B.H. Paw, A.M. Falick, H.E. Witkowska, J. Flint, D. Higgs, et al. 2003. Characterization of embryonic globin genes of the zebrafish. *Dev. Biol.* 255:48-61.
12. Glenn, N.O., J.A. Schumacher, H.J. Kim, E.J. Zhao, J. Skerniskyte, and S. Sumanas. 2014. Distinct regulation of the anterior and posterior myeloperoxidase expression by ETV2 and GATA1 during primitive Granulopoiesis in zebrafish. *Dev. Biol.* 393:149-159.

Received 22 August 2014; accepted 16 September 2014.

Address correspondence to Saulius Sumanas, Division of Developmental Biology, Cincinnati Children's Hospital Medical Center, Cincinnati, OH. E-mail: saulius.sumanas@cchmc.org

To purchase reprints of this article, contact: biotechniques@fosterprinting.com

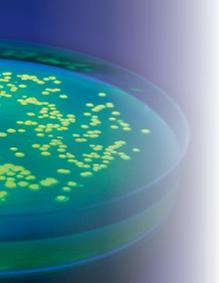

Plasmid & Minicircle DNA Service

- Customized plasmids and *minicircles*
- *InStock* reporter genes
- AAV helper & packaging plasmids

PlasmidFactory.com



The Minicircle Company



PlasmidFactory GmbH & Co. KG | Meisenstraße 96 | D-33607 Bielefeld | Germany | Fon ++49 521 2997350